

Determination of Pharmaceuticals by Capillary HPLC-MS/MS

Authors: S.V. Calugaru, E.P. Gates, W.R. West, M.L. Lee

Abstract

This application note demonstrates the use of a compact portable capillary liquid chromatograph, the Axcend Focus LC, coupled to an Agilent Ultivo triple quadrupole mass spectrometer for quantitative analysis of pharmaceutical drugs in model aqueous samples. Two test mixtures were analyzed in this study. Mix #1 contained acetaminophen, caffeine, carbamazepine, ciprofloxacin, erythromycin, fluoxetine, sulfamethoxazole, and trimethoprim. Mix #2 was composed of gemfibrozil, ibuprofen, naproxen, and triclosan. All compounds in the test mixtures were identified and calibration curves were constructed in the concentration ranges of 30-1,000 ng/mL for carbamazepine, 50-10,000 ng/mL for trimethoprim, 100-10,000 ng/mL for acetaminophen, 100-5,000 ng/mL for caffeine, 300-10,000 ng/mL for fluoxetine, 300-3,000 ng/mL for sulfamethoxazole, 500-10,000 ng/mL for erythromycin and ciprofloxacin, 30-10,000 ng/mL for gemfibrozil, ibuprofen and naproxen, and 30-1,000 ng/mL for triclosan.

Introduction

The widespread use of pharmaceuticals including hormones, antibiotics, non-steroidal anti-inflammatory drugs, antidepressants, and antifungal agents results in significant discharge of these drugs and their metabolites into wastewater [1]. Concerns have been raised over their potential to contaminate ground and surface waters, and to enter drinking water supplies. Therefore, monitoring pharmaceuticals in aqueous samples represents an important analytical task. Using capillary LC for this purpose provides many advantages including dramatically reduced consumption of toxic and expensive organic solvents [2] and potentially higher sensitivity when using MS detection resulting from intrinsically low mobile phase flow rates [3].

Experimental Details

Instrumentation

A compact capillary LC, the Axcend Focus LC with 2.2.0 Axcend Drive software (Axcend, Provo, Utah, USA), was used in this study. The LC was interfaced to a triple quadrupole mass spectrometer, the Ultivo LC/TQ G6465B equipped with Jet Stream Electrospray Ionization source (Agilent Technologies, Santa Clara, California, USA). MassHunter software (Acquisition: v1.1, Qualitative Analysis: v10.0, Quantitative Analysis: v10.0) was used for MS control and data processing. To accommodate the capillary LC mobile phase microflow rates, the regular nebulizer of the ion source was replaced with a microflow nebulizer (Part No G1946-67260) obtained from Agilent. A 25 cm long, 360 μm OD, 25 μm ID PEEKsil tubing (Part No. 0624374, Trajan, Melbourne, Victoria, Australia) served as a transfer line. One end of this tubing was attached to the solvent outlet from the column cartridge using a PEEK zero-dead-volume 360 μm union with a 50 μm bore hole (Part No. C360UPK2, VICI Valco Instruments, Houston, Texas, USA). The other end of the tubing was connected to the microflow nebulizer using a stainless steel 1/16" to 360 μm zero-dead-volume reducing union with a 100 μm bore hole (Part No. C360RUS64, VICI Valco Instruments). Direct sample infusion into the mass spectrometer for optimization of acquisition parameters was performed with a Model 22 syringe pump from Harvard Apparatus (Holliston, Massachusetts, USA).

APPLICATION NOTE

Determination of Pharmaceuticals by Capillary HPLC-MS/MS

Chemicals and Solvents

Mixtures of pharmaceuticals (Mix #1 containing acetaminophen, caffeine, carbamazepine, ciprofloxacin, erythromycin, fluoxetine, sulfamethoxazole, and trimethoprim at 200 µg/mL each in methanol, and Mix #2 containing gemfibrozil, ibuprofen, naproxen, and triclosan at 200 µg/mL each in methanol) were purchased from Restek (Bellefonte, Pennsylvania, USA). Calibration standards were prepared by serial dilutions of these test mixtures in water. D4-Acetaminophen (100 µg/mL in methanol) and D10-carbamazepine (100 µg/mL in methanol) were purchased from HPC Standards (Atlanta, Georgia, USA) and used as internal standards in Mix #1 samples. D3-Ibuprofen (100 µg/mL in methanol) was purchased from Cerilliant (Round Rock, Texas, USA) and used as an internal standard in Mix #2 samples. LC-MS grade solvents (water and acetonitrile) were obtained from Sigma-Aldrich (MilliporeSigma, St. Louis, Missouri, USA) and LC-MS grade formic acid was purchased from ThermoFisher Scientific (Waltham, Massachusetts, USA).

HPLC Methods

Mix #1 was analyzed using an HSS T3 Acquity UPLC column (10 cm x 150 µm i.d., 1.8 µm particle size) from Waters (Milford, Massachusetts, USA). A small diameter flow-cell detector was used for monitoring UV absorption at 275 nm. Mix #2 was analyzed using a C18 column (10 cm x 150 µm i.d., 1.7 µm particle size) from CoAnn Technologies (Richland, Washington, USA). UV absorption was monitored at 235 nm using an on-column detector. In both LC methods, solvent A was composed of 97% water, 3% acetonitrile, and 0.1% formic acid, and solvent B was composed of 97% acetonitrile, 3% water, and 0.1% formic acid. The flow rate was 1 µL/min and the injection volume was 250 nL (full loop). For the separation of Mix #1 components, the content of solvent B in the mobile phase was 3% during the initial 0.5 min after sample injection, then it was linearly ramped to 26% over 1.5 min, held at 26% for 5 min, linearly ramped to 97% over the next 3 min, and finally it was held at 97% for 5 min. For Mix #2 the content of solvent B in the linear mobile phase gradient was 3% for the initial 0.5 min, then it was increased to 60% over 1 min, to 77% over the next 2.5 min, to 97% over the following 1.5 min, and finally it was held at 97% for 3 min.

MS Methods

Multiple reaction monitoring (MRM) was used in positive and negative polarity modes for the analysis of both test mixtures. Deuterated acetaminophen, carbamazepine and Ibuprofen were used as internal standards. The MRM methods were optimized using MassHunter optimizer software by injecting the respective diluted test mixture (each component at 100 ng/mL) into the ion source at a flow rate of 1 µL/min with a syringe pump. The identified optimum fragmentor voltage and collision energy values for the quantifier and qualifier ions of the individual analytes and internal standards are presented in Table 1 for Mix #1 and in Table 2 for Mix #2. The other MS parameters were as follows: capillary voltage, 3000 V; gas temperature, 200 °C; gas flow rate, 5 L/min; nebulizer pressure, 10 psi; and dwell time, 50 ms except for fluoxetine and sulfamethoxazole, which was 100 ms.

Results and Discussion

Identification of Analytes

Two transitions between the parent ion and the product ions were selected for monitoring analytes and internal standards. The first transition was used for analyte quantification and the second one for qualification (Tables 1 and 2). No fragmentation of triclosan was observed under the tested conditions. For this pharmaceutical, only the molecular ion was monitored (Table 2). The chemical structures of the analyzed pharmaceuticals are shown in Figure 1 for Mix #1 and Figure 2 for Mix #2. The chromatographic separations of the components of Mix #1 and Mix #2 with UV detection at 275 nm and 235 nm are presented in Figures 3 and 4, respectively. Each analyte in the UV chromatograms was identified based on MS data. The extracted ion chromatograms for Mix #1 and Mix #2 are shown in Figures 5 and 6, respectively.

APPLICATION NOTE

Determination of Pharmaceuticals by Capillary HPLC-MS/MS

Quantification of Analytes

To normalize the intensities of the MS signals from Mix #1 components, deuterated analogs of acetaminophen (D4-acetaminophen) and carbamazepine (D10-carbamazepine) were used as internal standards at a final concentration of 500 ng/mL each. D4-acetaminophen served as an internal standard for acetaminophen and ciprofloxacin, while D10-carbamazepine was used for the other compounds in Mix #1. Linear calibration curves with a regression coefficient $R^2 > 0.99$ were obtained in the following concentration ranges: 30-1,000 ng/mL for carbamazepine, 50-10,000 ng/mL for trimethoprim, 100-10,000 ng/mL for acetaminophen, 100-5,000 ng/mL for caffeine, 300-10,000 ng/mL for fluoxetine, 300-3,000 ng/mL for sulfamethoxazole, and 500-10,000 ng/mL for erythromycin. Quadratic regression was used to obtain calibration curves with a regression coefficient $R^2 > 0.99$ in the concentration range of 500-10,000 ng/mL for ciprofloxacin. Typical calibration curves are shown in Figures 7 and 8 for acetaminophen and carbamazepine, respectively. In the case of Mix #2, D3-ibuprofen was used at a final concentration of 250 ng/mL as an internal standard for all analytes. Calibration curves with a regression coefficient $R^2 > 0.99$ were obtained in the following concentration ranges: 30-10,000 ng/mL for gemfibrozil, ibuprofen, and naproxen (linear regression), and 30-1,000 ng/mL for triclosan (quadratic regression). Representative calibration curves are shown in Figures 9 and 10 for ibuprofen and naproxen, respectively.

Table 1. Optimized Acquisition Parameters for MS Detection of Precursor and Product Ions of Compounds in Pharmaceutical Mix #1

Compound Name	Precursor Ion (m/z)	Quantifier Product Ion (m/z)	Qualifier Product Ion (m/z)	Fragmentor Voltage (V)	Collision Energy (eV)	Polarity
Acetaminophen	150.1	106.9	-	94	25	Negative
Acetaminophen	152.1	-	109.8	100	25	Positive
Caffeine	195.1	138.1	109.8	98	21	Positive
Ciprofloxacin	332.1	230.8	314.4	100	35	Positive
D4-Acetaminophen	154.1	110.9	-	98	17	Negative
D4-Acetaminophen	156.1	-	114.1	75	10	Positive
Trimethoprim	291.2	80.8	122.9	100	50	Positive
Carbamazepine	237.1	194.1	179.1	100	25	Positive
D10-Carbamazepine	247.2	204.2	186	105	25	Positive
Erythromycin	734.5	576.5	158	100	25	Positive
Fluoxetine	310.1	91.3	116.9	100	55	Positive
Sulfamethoxazole	254.1	64.9	91.8	100	45	Positive

Table 2. Optimized Acquisition Parameters for MS Detection of Precursor and Product Ions of Compounds in Pharmaceutical Mix #2

Compound Name	Precursor Ion (<i>m/z</i>)	Quantifier Product Ion (<i>m/z</i>)	Qualifier Product Ion (<i>m/z</i>)	Fragmentor Voltage (V)	Collision Energy (eV)	Polarity
Gemfibrozil	251.1	83.1	-	50	10	Positive
Gemfibrozil	249.0	-	120.9	100	0	Negative
Ibuprofen	205.0	160.9	-	100	0	Negative
Ibuprofen	207.1	-	118.9	50	10	Positive
D3-Ibuprofen	210.1	164.1	-	100	5	Positive
D3-Ibuprofen	208.1	163.9	-	100	5	Negative
Naproxen	229.0	184.9	-	100	0	Negative
Naproxen	231.1	-	185.1	50	10	Positive
Triclosan	286.9	286.9	-	100	0	Negative

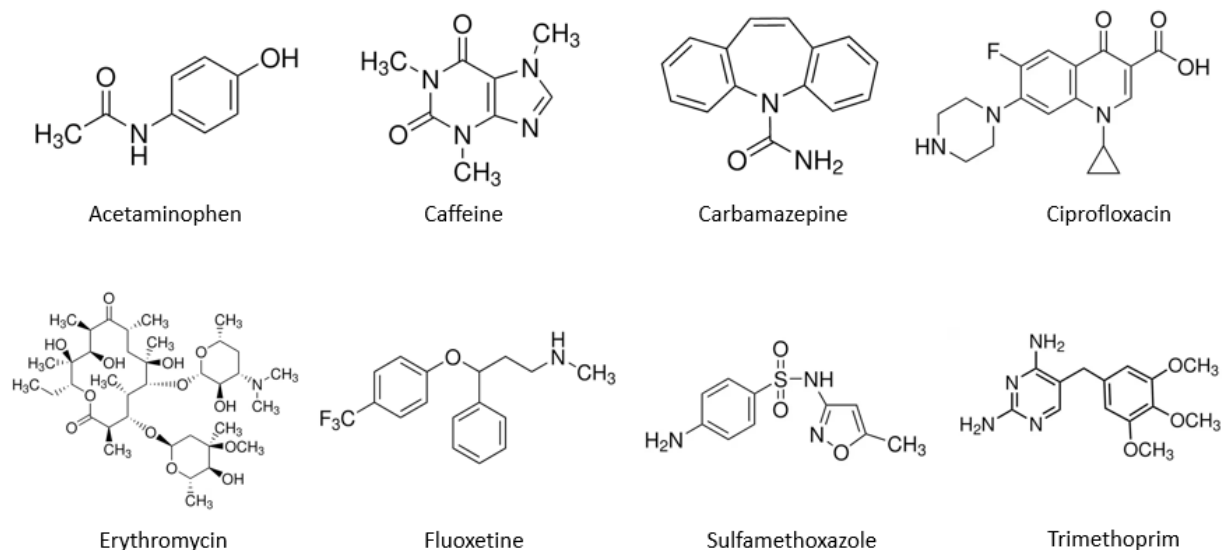


Figure 1. Chemical structures of the compounds in Mix #1.

APPLICATION NOTE
 Determination of Pharmaceuticals
 by Capillary HPLC-MS/MS

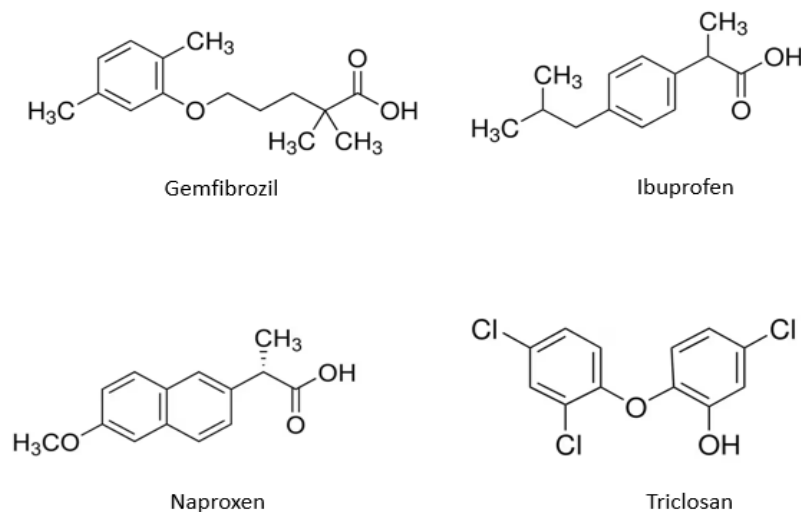


Figure 2. Chemical structures of the compounds in Mix #2.

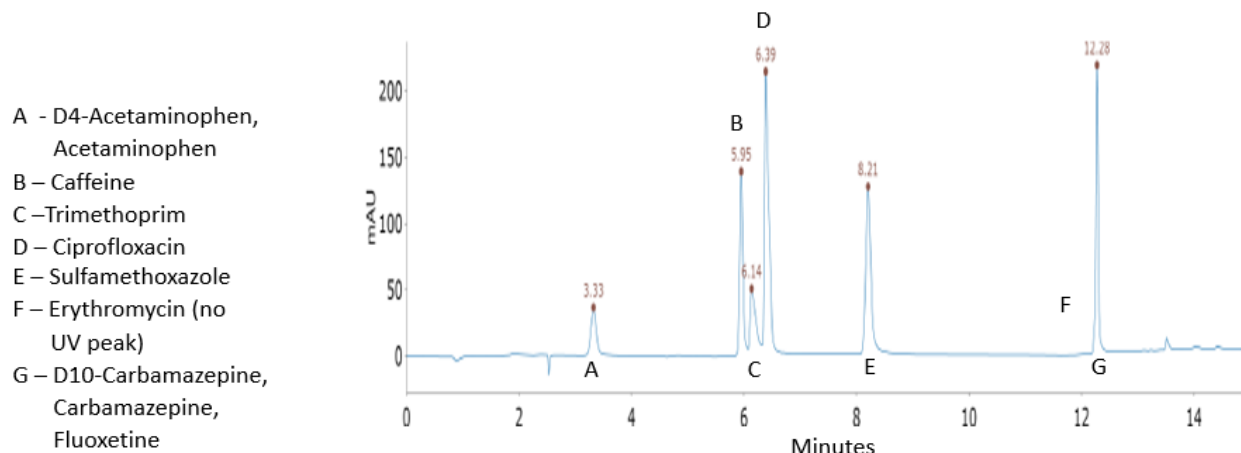


Figure 3. Analysis of Mix #1 (5 µg/mL each analyte) using UV detection at 275 nm. Identifications are based on MS data (see Figure 5).

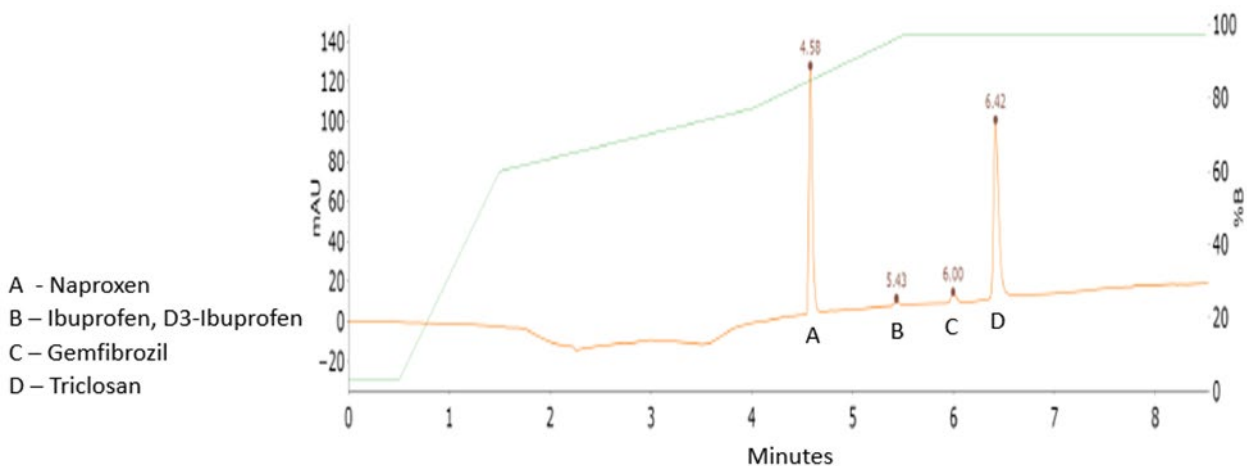


Figure 4. Analysis of Mix #2 (10 µg/mL each analyte) using UV detection at 275 nm. Identifications are based on MS data (see Figure 6).

APPLICATION NOTE
Determination of Pharmaceuticals
by Capillary HPLC-MS/MS

- A – D4-Acetaminophen,
Acetaminophen
- B – Caffeine
- C – Trimethoprim
- D – Ciprofloxacin
- E – Sulfamethoxazole
- F – Erythromycin (no UV peak)
- G – D10-Carbamazepine,
Carbamazepine,
Fluoxetine

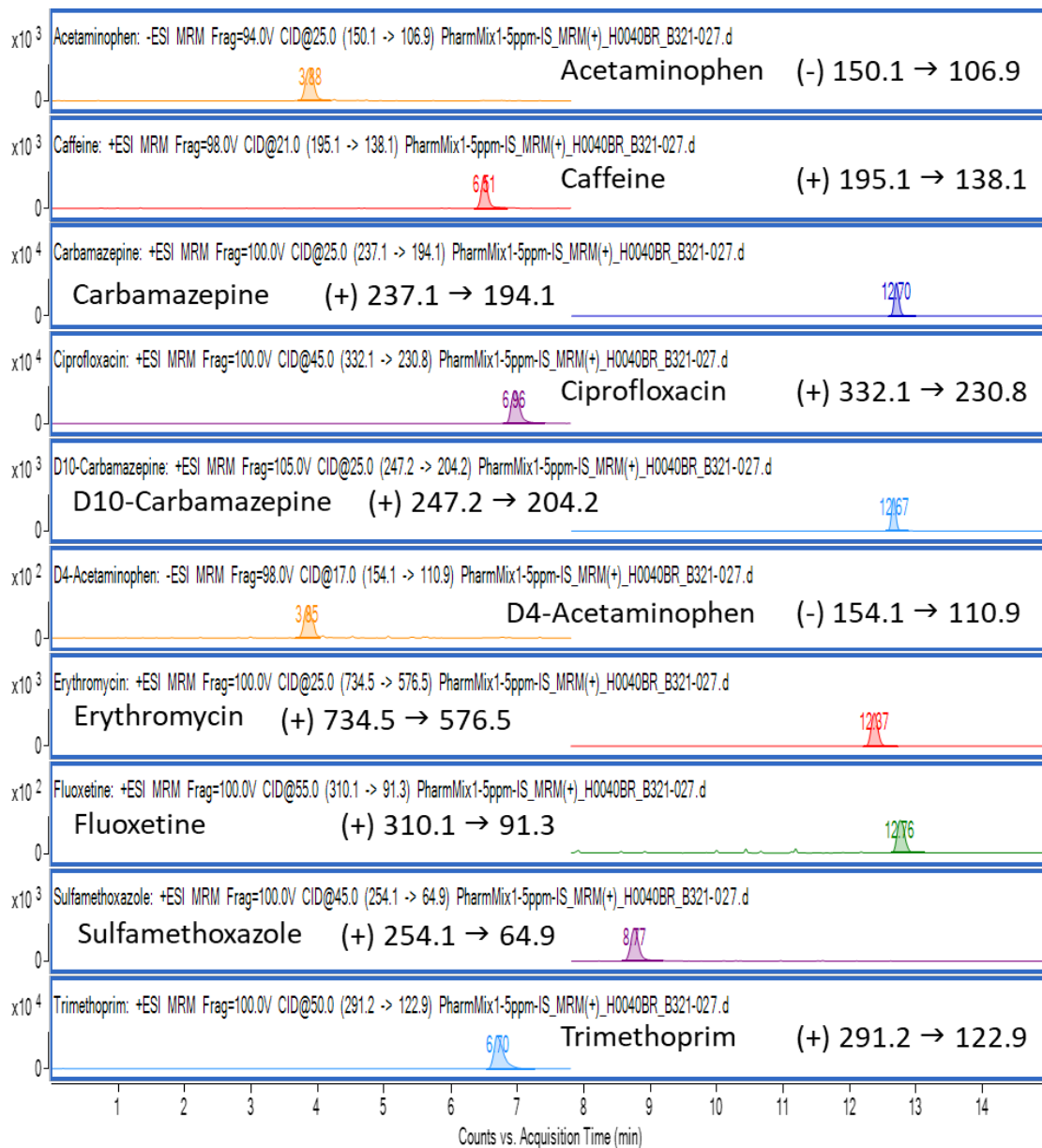
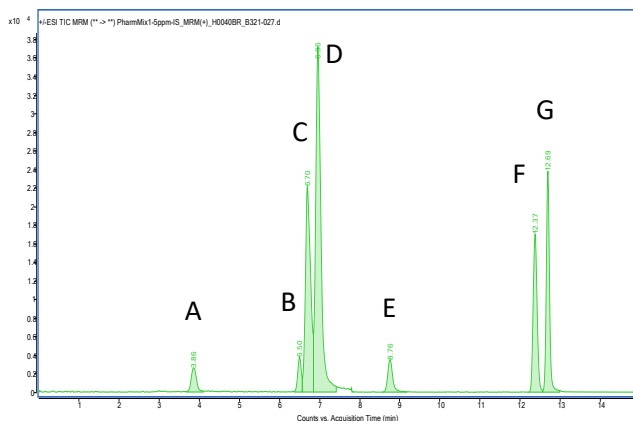


Figure 5. Analysis of Mix #1 (5 µg/mL each analyte) using MS detection in the MRM mode: total ion current (top) and extracted ion chromatograms (bottom).

APPLICATION NOTE
Determination of Pharmaceuticals
by Capillary HPLC-MS/MS

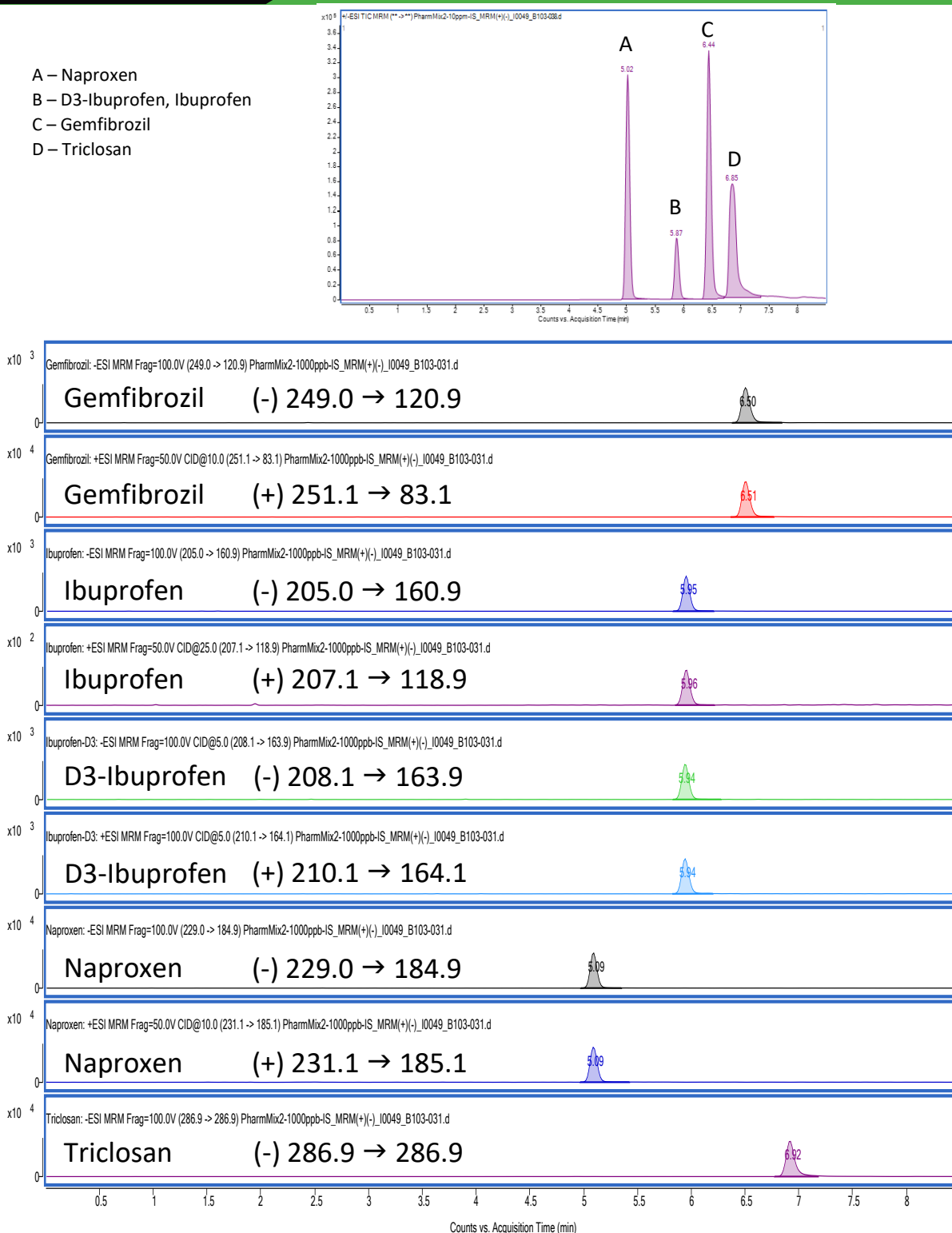


Figure 6. Analysis of Mix #2 (10 µg/mL each analyte) using MS detection in the MRM mode: total ion current (top) and extracted ion chromatograms (bottom).

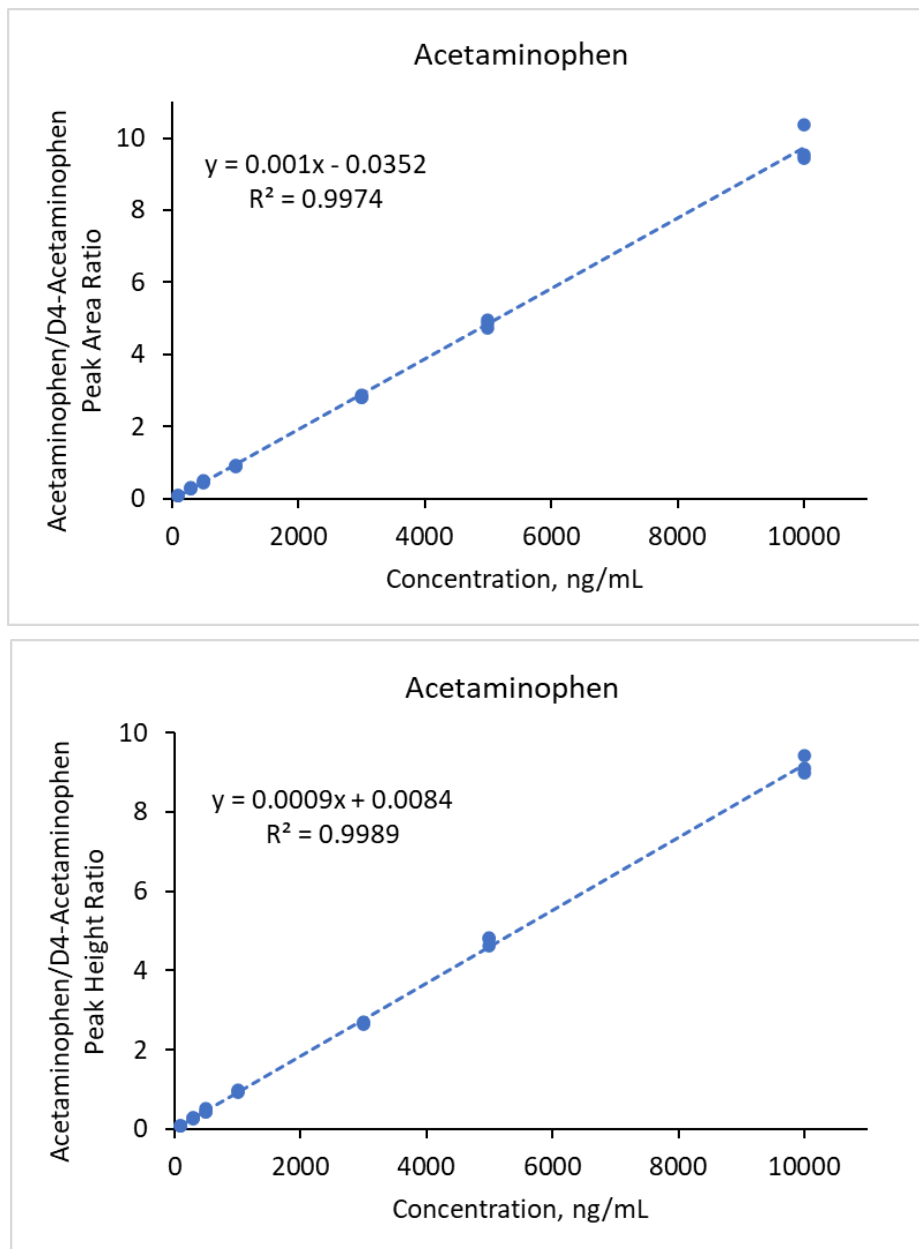


Figure 7. Calibration curves for acetaminophen in the concentration range of 100-10,000 ng/mL based on peak areas (top) and peak heights (bottom).

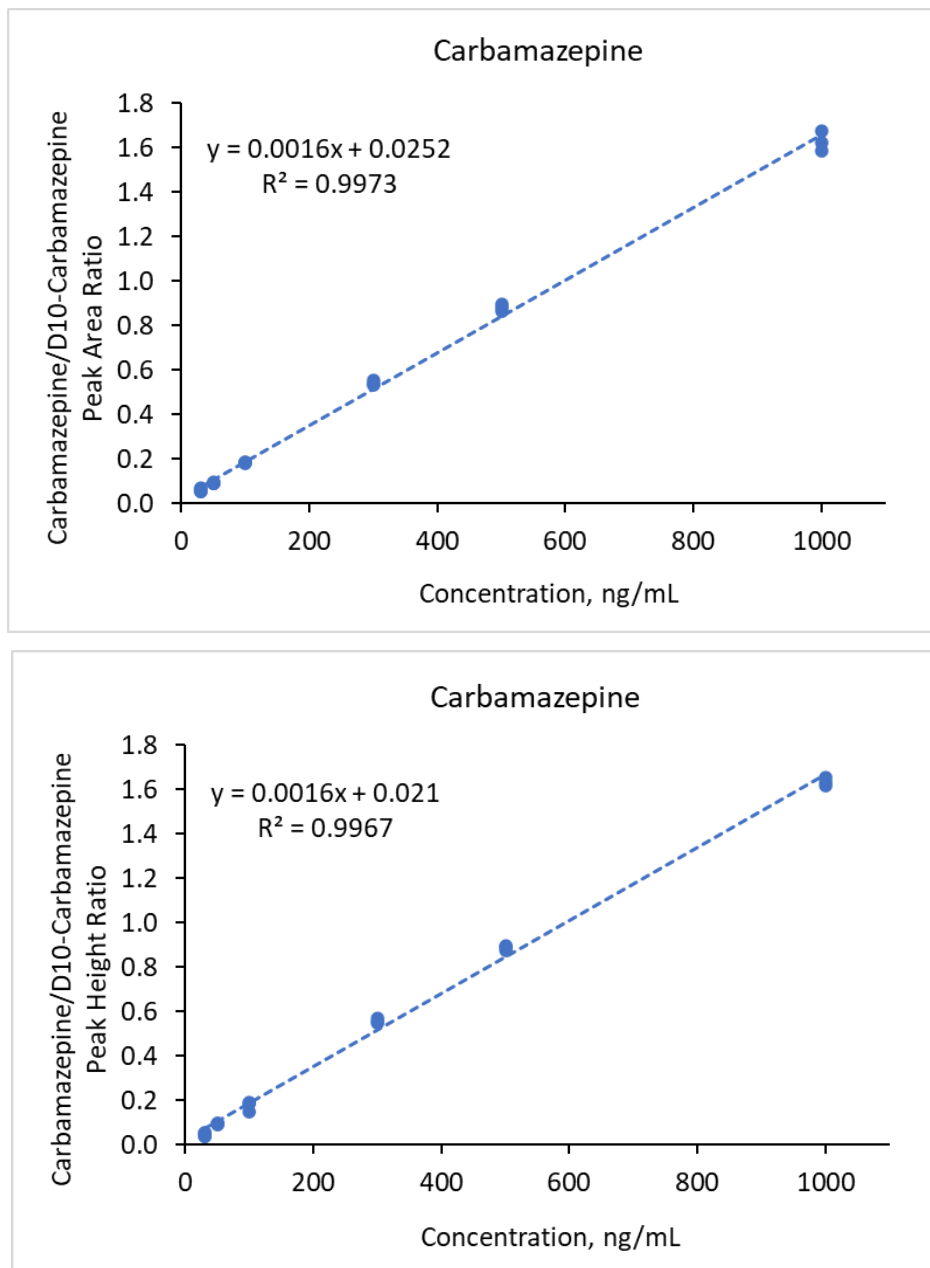


Figure 8. Calibration curves for carbamazepine in the concentration range of 30-1,000 ng/mL based on peak areas (top) and peak heights (bottom).

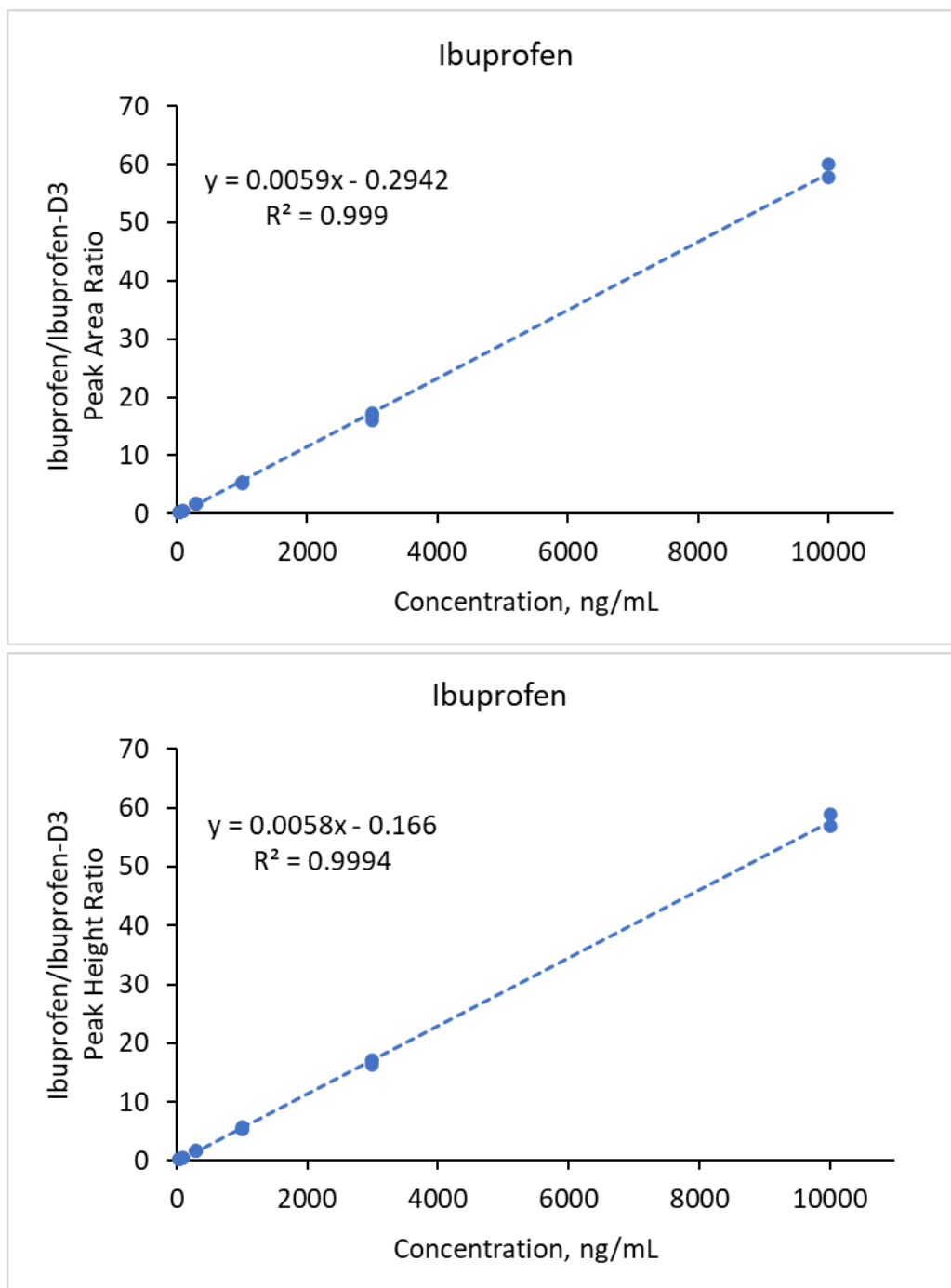


Figure 9. Calibration curves for ibuprofen in the concentration range of 30-10,000 ng/mL based on peak areas (top) and peak heights (bottom).

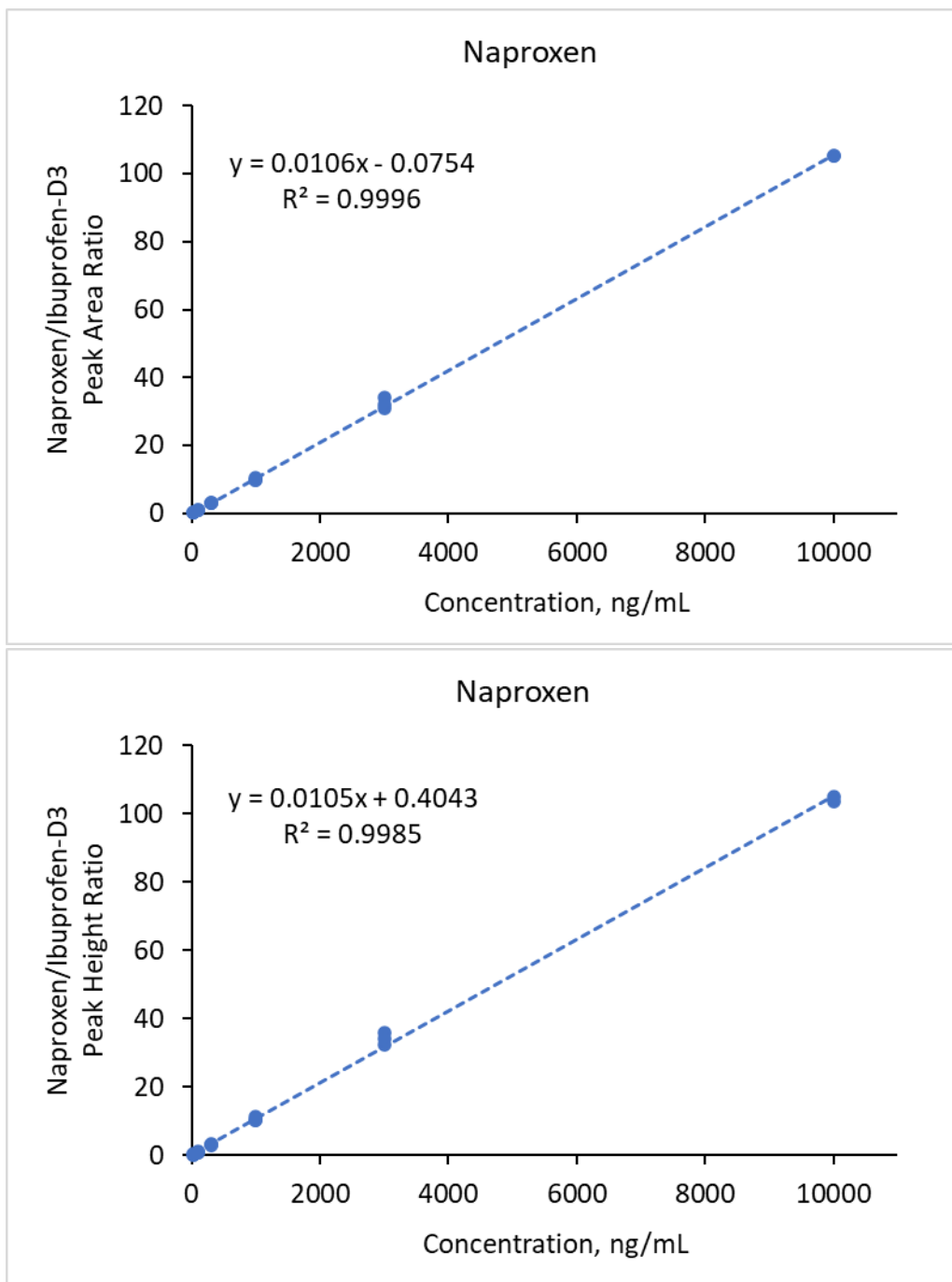


Figure 10. Calibration curves for naproxen in the concentration range of 30-10000 ng/mL based on peak areas (top) and peak heights (bottom).

APPLICATION NOTE

Determination of Pharmaceuticals by Capillary HPLC-MS/MS

The results of this study can be used in the development of quantitative methods for the analysis of pharmaceuticals in complex matrices. The presented HPLC-MS/MS methods can be applied, for example, to the analysis of various fractions in the development and optimization of sample preparation procedures.

Conclusions

Identification and quantification of pharmaceuticals in mixtures of components in model aqueous samples based on capillary HPLC coupled to tandem MS was successfully performed.

Calibration curves with a regression coefficient R^2 above 0.99 were obtained in the concentration ranges typically covering over two orders of magnitude (for example, 30-10,000 ng/mL for acetaminophen).

The developed HPLC-MS/MS methods can serve as a basis for further development of quantitative capillary HPLC-MS/MS determination of pharmaceuticals in biological matrices.

References

1. M.R. Letsoalo, T. Sithole, S. Mufamadi, Z. Mazhandu, M. Sillanpaa, A. Kaushik, T. Mashifana. Efficient Detection and Treatment of Pharmaceutical Contaminants to Produce Clean Water for Better Health and Environment. *J Cleaner Production*, 2023, 387, 135798.
2. S.W. Foster, X. Xie, M. Pham, P.A. Peaden, L.M. Patil, L.T. Tolley, P.B. Farnsworth, H.D. Tolley, M.L. Lee, J.P. Grinias. Portable Capillary Liquid Chromatography for Pharmaceutical and Illicit Drug Analysis. *J Sep Sci*, 2020, 43 (9-10), 1623–1627.
3. D.A.V. Medina, E.V.S. Maciel, A.L. de Toffoli, F.M. Lanças. Miniaturization of Liquid Chromatography Coupled to Mass Spectrometry 2. Achievements on Modern Instrumentation for Miniaturized Liquid Chromatography Coupled to Mass Spectrometry. *TrAC Trends in Analytical Chemistry*, 2020, 128, 115910.



Axcend
5252 N. Edgewood Dr. #185
Provo, UT 84604 USA
801-405-9545

www.axcendcorp.com

Axcend, the Axcend logo, and Axcend Focus LC are each trademarks of Axcend. All other trademarks are property of their respective owners. Copyright 2024 © Axcend